

A new twist on the nail:  
revolutionizing dynamic  
compression with disruptive  
NiTiNOL technology




# DYNANAIL HELIX™

ACTIVE, ADAPTIVE HEALING FOR FOOT & ANKLE FUSION

## NAIL SYSTEM

### FEATURES & BENEFITS:

- Internal NiTiNOL element offers between 1 - 5 mm of active compression in response to bone settling or resorption.
- Available in headed and headless options for procedural and anatomical versatility. Sizes range from 50-100 mm in 5 mm increments.
- Threads allow for easy insertion, similar to a screw.
- NiTiNOL Compressive Element comes pre-stretched and pre-loaded on a disposable Nail Guide for enhanced efficiency.

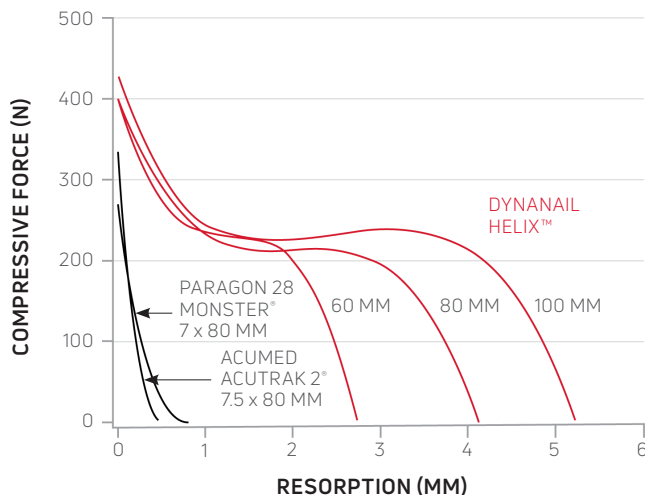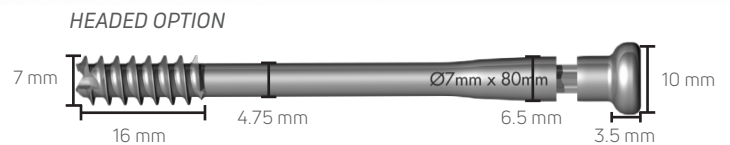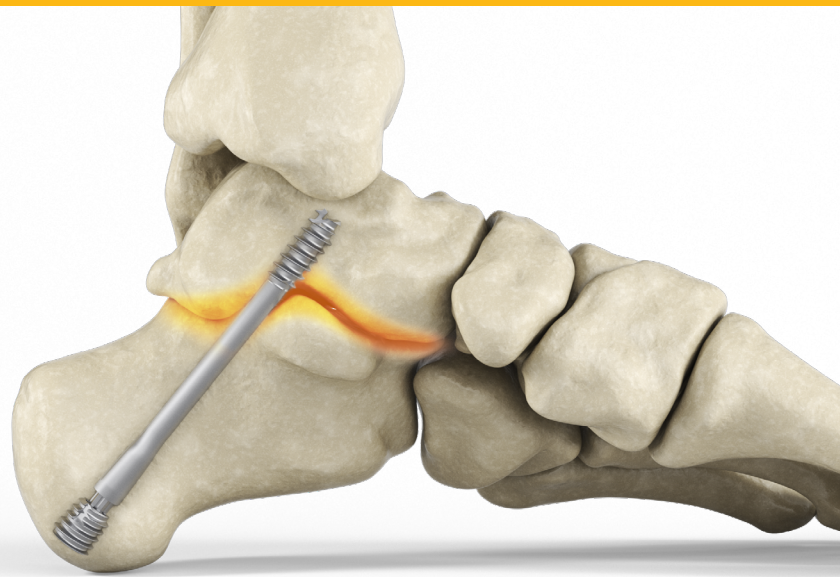

### SUSTAINED DYNAMIC COMPRESSION

Next generation technology encourages active, adaptive healing through sustained dynamic compression.

Compared to traditional static screws, DynaNail Helix provides significantly more compression in response to bone settling or resorption.<sup>1</sup>

Data represents loss of compression as function of increased bone resorption in synthetic bone. Laboratory testing does not necessarily indicate clinical performance.


Monster® is a registered trademark of Paragon 28, Inc

Acutrak 2® is a registered trademark of Acumed, LLC

## INSTRUMENTATION HIGHLIGHTS


### PRE-LOADED ASSEMBLY

DynaNail Helix™ is provided with the NiTiNOL Element pre-stretched and pre-loaded on a disposable Nail Guide that connects directly to the ratcheting handle for easy insertion.

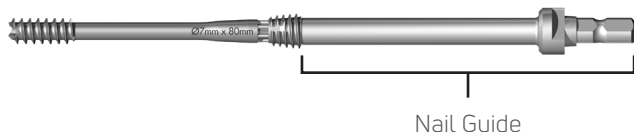

### PARALLEL PIN GUIDE

The Parallel Pin guide can be used to assist with parallel pin placement when using multiple implants.

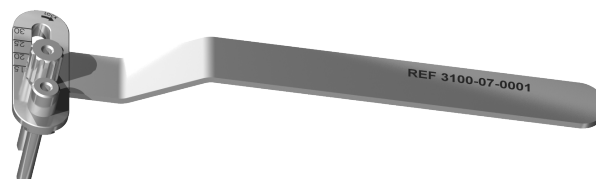


### DYNANAIL HELIX™ IMPLANTS

PART #	DESCRIPTION
3200-00-7050	Ø7 MM x 50 MM, HEADED FUSION NAIL
3200-00-7055	Ø7 MM x 55 MM, HEADED FUSION NAIL
3200-00-7060	Ø7 MM x 60 MM, HEADED FUSION NAIL
3200-00-7065	Ø7 MM x 65 MM, HEADED FUSION NAIL
3200-00-7070	Ø7 MM x 70 MM, HEADED FUSION NAIL
3200-00-7075	Ø7 MM x 75 MM, HEADED FUSION NAIL
3200-00-7080	Ø7 MM x 80 MM, HEADED FUSION NAIL
3200-00-7085	Ø7 MM x 85 MM, HEADED FUSION NAIL
3200-00-7090	Ø7 MM x 90 MM, HEADED FUSION NAIL
3200-00-7095	Ø7 MM x 95 MM, HEADED FUSION NAIL
3200-00-7100	Ø7 MM x 100 MM, HEADED FUSION NAIL
3200-01-7050	Ø7 MM x 50 MM, HEADLESS FUSION NAIL
3200-01-7055	Ø7 MM x 55 MM, HEADLESS FUSION NAIL
3200-01-7060	Ø7 MM x 60 MM, HEADLESS FUSION NAIL
3200-01-7065	Ø7 MM x 65 MM, HEADLESS FUSION NAIL
3200-01-7070	Ø7 MM x 70 MM, HEADLESS FUSION NAIL
3200-01-7075	Ø7 MM x 75 MM, HEADLESS FUSION NAIL
3200-01-7080	Ø7 MM x 80 MM, HEADLESS FUSION NAIL
3200-01-7085	Ø7 MM x 85 MM, HEADLESS FUSION NAIL
3200-01-7090	Ø7 MM x 90 MM, HEADLESS FUSION NAIL
3200-01-7095	Ø7 MM x 95 MM, HEADLESS FUSION NAIL
3200-01-7100	Ø7 MM x 100 MM, HEADLESS FUSION NAIL

### DYNANAIL HELIX™ SINGLE-USE INSTRUMENTS

PART #	DESCRIPTION
2900-04-0229	GUIDEWIRE, 2.4 MM x 229 MM
3100-01-0005	5.0 MM DRILL
3100-01-0007	PROFILE DRILL, 6.5 MM
3100-04-0001	TRAJECTORY WIRE, 1.0 x 260 MM
2201-09-0025	2.5 MM DRILL

### CHECK OUT MORE OF THE DYNANAIL FAMILY

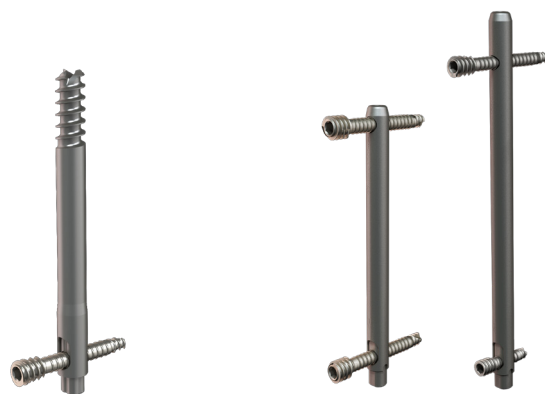

DYNANAIL HYBRID™  
FUSION SYSTEM

DYNANAIL MINI®  
FUSION SYSTEM

**enovis**

T 800.495.2919 F 877.778.3864

Medshape, Inc.  
1575 Northside Drive NW | Suite 440 | Atlanta, GA 30318 | U.S.A.  
enovis.com/foot-and-ankle

Copyright © 2022 DJO Foot and Ankle  
MK-10267 Rev 01

1. Data on File, DJO Global. 2022

Individual results may vary. Neither Medshape, Inc. nor any of the Enovis companies dispense medical advice. The contents of this document do not constitute medical, legal, or any other type of professional advice. Rather, please consult your healthcare professional for information on the courses of treatment, if any, which may be appropriate for you.